

TRADITIONAL EXAM ROOM WORKFLOW A



Some of the leading factors to consider in today's exam room design are improving patient care, accessibility, efficiency and workflow. This exam room provides the foundation for a basic caregiver-patient workflow, with separate zones for patient interviews and exams, and the option for additional privacy with public/private zones. While this traditional setup allows for immediate patient care, historical patient data and any additional medical data that supports patient care is not immediately available in the room, so caregivers must leave the room to retrieve this information.




TRADITIONAL EXAM ROOM WORKFLOW A




For better access to the patient, the examination table is positioned at an angle, providing access from three sides. The exam room is organized into four zones, including entry, patient care, family/visitor and organized storage with an integrated sink.

A lowered writing surface is adjacent to the storage zone, as is a chair for the initial patient interview. The overall room size is 10 ft. x 10 ft. with the entry door reverse-hinged for patient privacy.

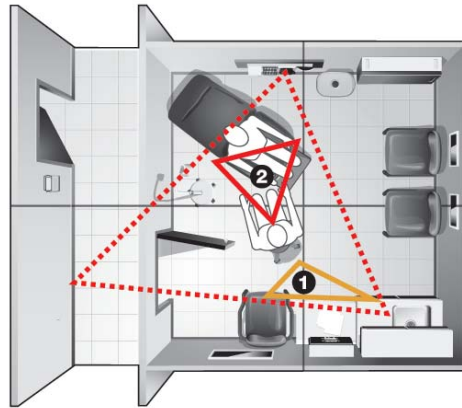

Traditional Workflow

Patient care is divided into separate zones for the patient interview **1** and patient exam area **2**. While the traditional exam room supports immediate patient care, many times historic patient data and additional medical data that support patient care are retrieved from outside the exam room.

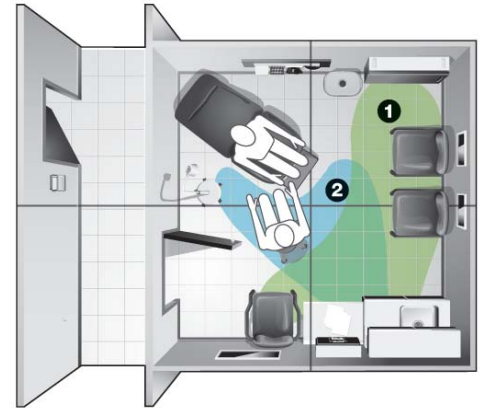

Public/Private Zones

This setup provides a public zone with side chairs to seat family/visitors, as well as a dressing nook for patients **1**. The patient care zone, also known as the private zone, allows access to equipment and supplies **2**. Space is provided for the initial patient interview, including a work surface to review paper records.

Midmark is an ISO 13485 and ISO 9001 Certified Company

For more information or a demonstration, contact your Midmark dealer or call:
1-800-MIDMARK
Fax: 1-800-365-8631
Outside the U.S.A. call: 1-937-526-3662
Fax: 1-937-526-8392
or visit our website at midmark.com

© 2012 Midmark Corporation
60 Vista Drive P.O. Box 286
Versailles, Ohio 45380-0286
Products subject to improvement
changes without notice
007-0840-00 (08/12)